

RollerBottle Standard digital

Description

Flexible roller mixer for bottles in various sizes. With robust roller drive without up and down.


Depending on the size of the bottles, individual rollers can be removed or moved (except the drive roller).

One roller is always the drive roller and is fixed. All other rollers can be moved by the user. RollerBottle offers space for a total of five rollers in the standard version.


Two rollers are installed in the standard version. Additional rollers are available as an option.

The force is transferred from the first to the other rollers by contact with the containers. The drive roller should be loaded with the largest and heaviest bottle.

The other rollers should be loaded with the same or less load. Therefore, please ensure that the contact from one roller to the next is not interrupted and that the containers are in contact.

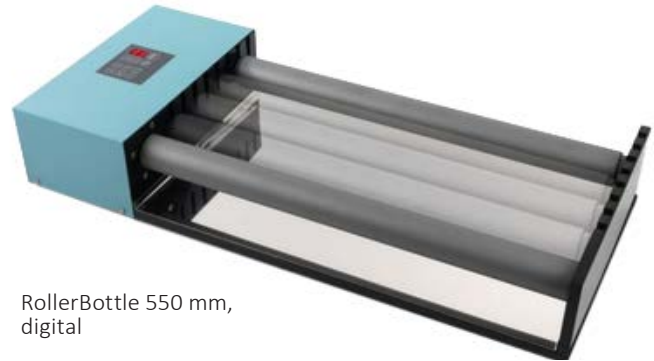

RollerBottle 550 mm, digital

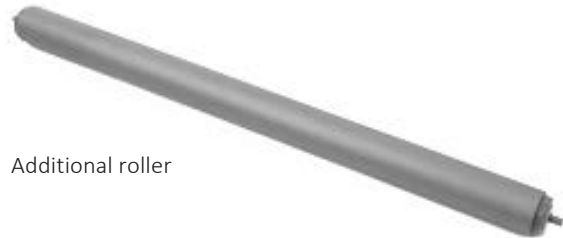

Additional roller

Optional Stacking Shelf

Stainless steel shelf suitable for RollerBottle. Sturdy shelf for stacking up to three appliances.




Example with three RollerBottles on stacking shelf


Model Standard digital

Length of rolls	Part No.
2 rolls à 550 mm	105.1011
2 rolls à 650 mm	105.1015

Other models with different roller lengths available on request.

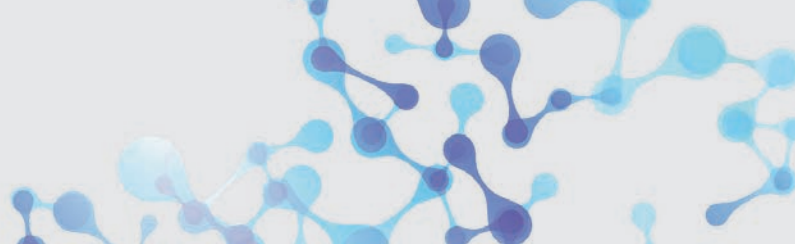

RollerBottle Standard digital

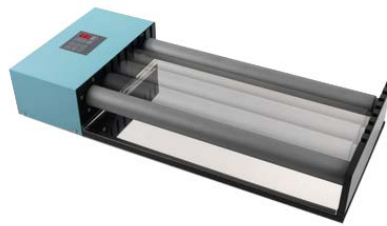

Technical Data

Number of rollers	max. 5
Max. loading	10 kg
Motor rating input	75 W
Motor rating output	50 W
Speed display	digital
Speed range	1 – 99 rpm
Operating time	1 sec – 99 h
Timer	Yes
Dimensions (W x D x H)	750 x 350 x 130 mm
Weight	10.5 kg
Standard equipment	2 rollers
∅ Roller	45.5 mm
Plate material	Steel with PE coated
Permissible ambient temperature	5 – 50 °C
Max. relative humidity	80 %
Voltage	230 V / 50/60 Hz

Roller Dimensions

The diameter of one roll is 45.5 mm.
If the rolls are placed directly next to each other,
there is 10.5 mm space between the individual rolls.
Measured from the center of the first roll to the
center of the second roll is 56 mm.