

IsoFreeze MCT Rack


Description

Ingenious optical control for your samples!
The middle section of the rack is filled with a thermochromic gel. This is able to change color depending on the temperature. This allows you to see immediately whether your samples are still at the right temperature.

IsoFreeze MCT racks ensure the protection of temperature-sensitive samples such as PCR reagents, cells or enzymes.

At 7°C there is a clear color change.
This allows you to see at a glance when the temperature of the samples rises above the original cooling temperature.
MCT racks keep the sample temperature below 4°C for up to four hours.

Each rack is supplied with an individual lid made of transparent polypropylene.

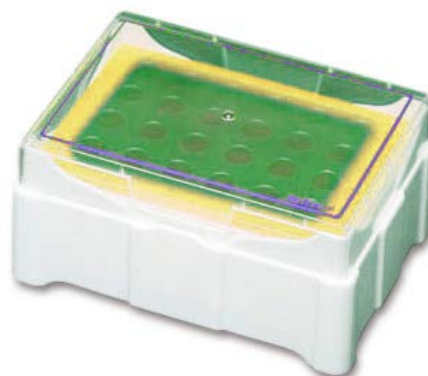

IsoFreeze MCT Rack green/yellow



IsoFreeze MCT Rack purple/pink

Technical Data

Capacity	24 Tubes with 0.5 / 1.5 / 2.0 ml volume
Dimensions (Rack + Lid)	146 x 105 x 70 mm (L x D x H)
Material	Polypropylene

Models

Description	Color	Packaging unit	Part number
IsoFreeze MCT Rack Suitable for tubes: 24 x 0,5 / 1,5 / 2,0 ml	green/yellow	1 rack with lid	104.0047
IsoFreeze MCT Rack Suitable for tubes: 24 x 0,5 / 1,5 / 2,0 ml	purple/pink	1 rack with lid	104.0044